

HEK-hERG Cell Line

Drug-induced blockage of the *human ether-a-go-go* (hERG) potassium channel is associated with a prolonged QT interval, which may predispose individuals to arrhythmias such as Torsades de Pointes. Therefore early screening of compounds for hERG channel activity is a vital step in the drug-discovery process.

Our HEK-hERG cell line was developed and validated for use as a screening tool for compounds which may lead to QT prolongation.

Our HEK-hERG cell line:

- **Stably expresses the hERG potassium channel**
- **Has been validated using electrophysiology and flux assay**
- **Is suitable for high-throughput screening (HTS)**

Electrophysiology

Electrophysiology experiments were conducted using standard patch clamp techniques.

Validation of HEK-hERG cell line was carried out using whole cell patch clamping to show the family of tail currents and I-V relationships (Figures 1-2).

As shown in Figure 3, the background HEK conductance is usually less than 1-200pA, maximal and there is no tail current. Under certain culture conditions one sees a shift in the voltage-dependence of activation of hERG such that significant current will be triggered with a prolonged depolarization at even -60mV.

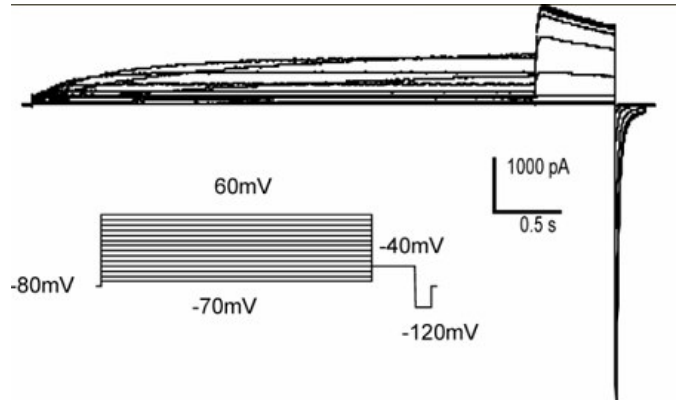


Fig. 1. hERG K⁺ channel activity transfected in HEK cells. Representative current tracings showing currents were elicited by depolarizing steps to various levels (between 40 and -80mV) followed by repolarizing steps to -70 mV.

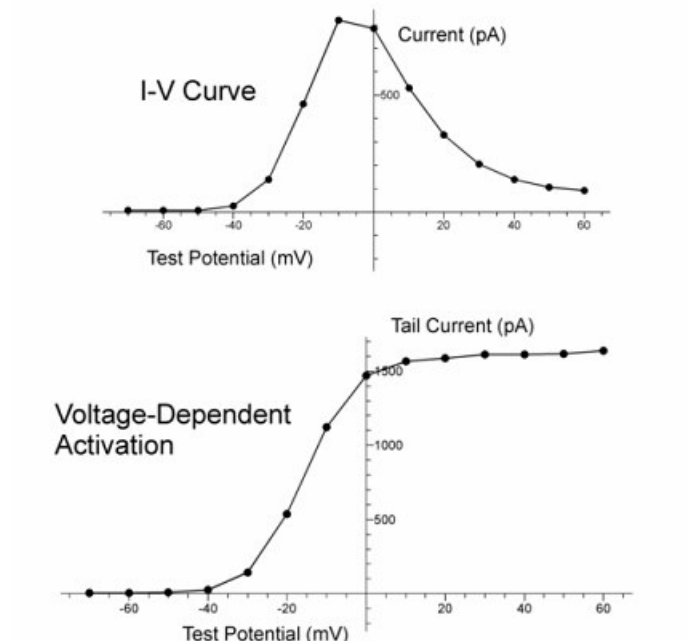


Fig. 2. hERG activity transfected in HEK cells showing I-V relationship, and voltage dependent activation. The voltage dependent activation curves plotted from peak tail currents during a repolarizing step to -70 mV after depolarizing to various voltages. Data are means and S.E., n = 8 per group.

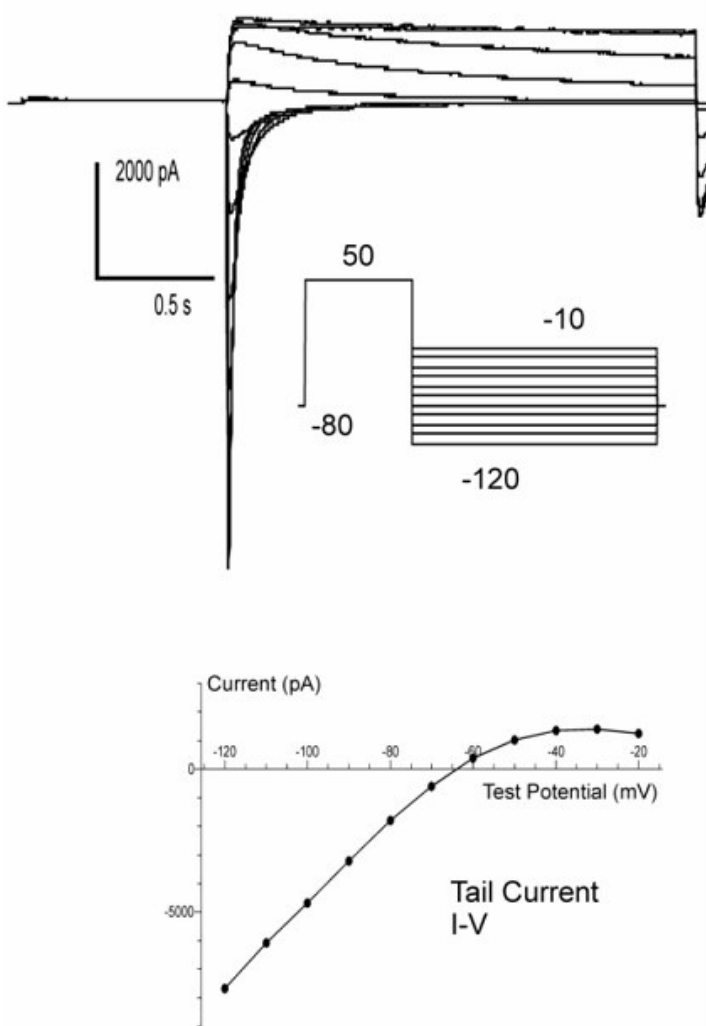


Fig. 3. Typical recording from a HEK-hERG cell using voltage-clamp protocols.

Rubidium Efflux Assay

The HEK-hERG cell line also validated for use with a cell based flux assay.

Using Aurora Biomed's Rubidium Efflux Assay protocol, the basal efflux was measured to be 5%. Following a 6 minute activation time using 60mM KCl, the maximal activation-induced efflux was 87, providing a more than 5-fold window of detection (Figure 4).

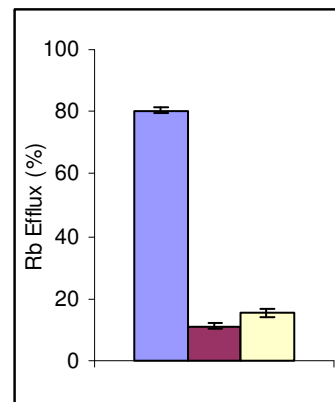


Fig. 4. Activation of hERG in the HEK cell line. An activation time of 6 minutes with 60 mM KCl leads to an 86% activation efflux (red bar), a 15% basal efflux (blue bar), and inhibition by 1 μM astemizole (yellow bar).

The concentration-response curves of three common hERG blockers, pimozone, astemizole and terfenadine, were determined using the HEK-hERG cell line. Figure 5 shows the IC₅₀ determinations of these three positive blockers. Table 1 shows the relative potency of these hERG blockers as determined using the HEK-hERG cell line and the Rubidium Efflux Assay.

Lastly, we calculated the concentration-response curves of normalized Rb⁺ flux data for these three commonly known hERG blockers, pimozone, astemizole and terfenadine, using a maximum concentration of 1 μM, 1 μM, and 3 μM respectively, as shown in Figure 6.

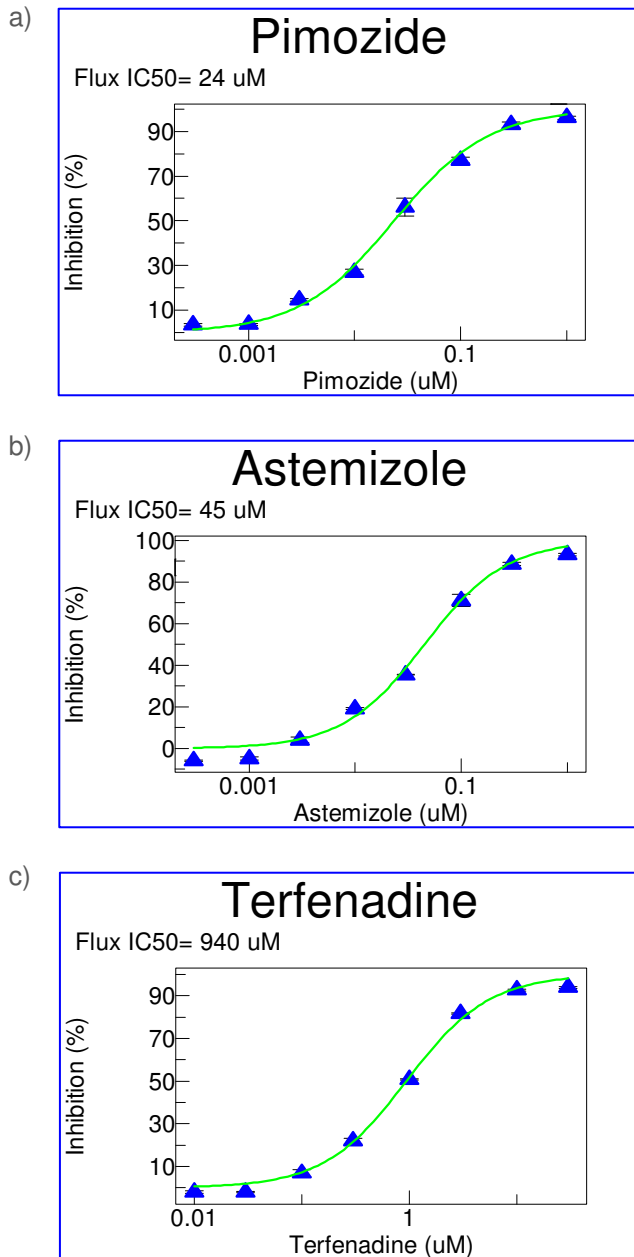


Fig. 5. IC₅₀ of (a) pimozide, (b) astemizole, and (c) terfenadine on HEK-hERG cell line as determined by the Rubidium Efflux Assay.

Table 1. Relative potencies of hERG blockers as determined by the Rubidium Flux Assay.

Drug	IC ₅₀ (μ M)
Pimozide	0.024
Astemizole	0.045
Terfenadine	0.940

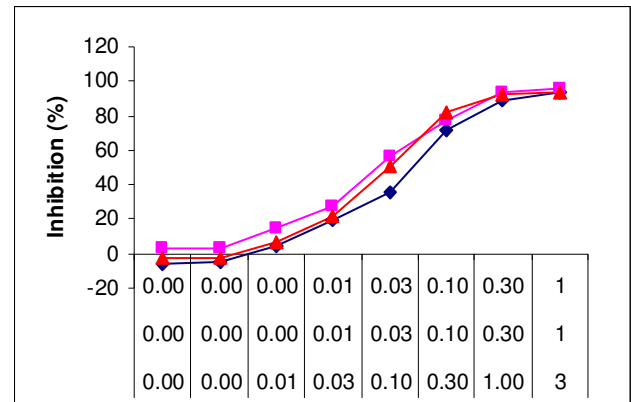


Fig. 6. The concentration-response curves of normalized Rb⁺ flux data of pimozide, astemizole, and terfenadine.